DIGITAL TERMOMETER BRUKSANVISNING

VARNING:

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DIGITAL THERMOMETER OWNER'S MANUAL

WARNING

- ⚠ Read the instructions carefully before using the digital thermometer.
- ⚠ Choking hazard: Thermometer lid and battery may be fatal if swallowed.
Keep it away from children and pets without adult supervision.
- ⚠ Intended use is for oral, rectal and armpit.
- ⚠ Do not place the thermometer battery near extreme heat.
- ⚠ It is recommended that performance be checked every two years.
- ⚠ Remove the battery from the unit when not in use for a long time.
- ⚠ Do not attempt to measure when the thermometer is wet, as incorrect readings may occur.
- ⚠ Do not bite the thermometer. This can lead to breakage and / or injury
- ⚠ Do not attempt to disassemble or repair the thermometer. This may result in inaccurate readings.
- ⚠ After each use, disinfect the thermometer.
- ⚠ Do not force the thermometer into the rectum. Stop insertion and interrupt the measurement if any pain occurs. Failure to do so may result in personal injury.
- ⚠ Do not use the thermometer orally after using it rectally.


PLEASE READ CAREFULLY BEFORE USE

This digital thermometer provides a quick and highly accurate reading of a individual's body temperature. Please read all instructions before use.

CONTENTS

1 Thermometer, 1 Owner's manual, 1 Storage Case

Product illustration

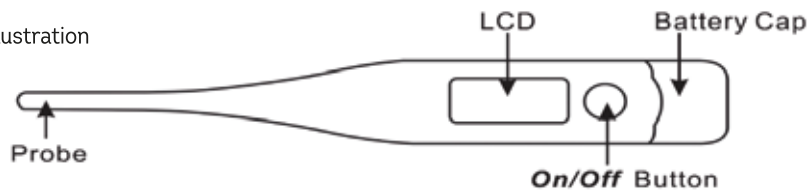PRECAUTION

The performance of the device may be degraded if one or more of the following occur:

- ⚠ Operation outside the manufacturer's stated temperature and humidity range.
- ⚠ Storage outside the manufacturer's stated temperature and humidity range.
- ⚠ Mechanical shock (For example, drop).
- ⚠ Patient temperature is below ambient temperature.

SYMBOL EXPLANATION

- | | | | |
|---|--|-----|--------------|
| → | Direct Current | LOT | Batch code |
| | Type BF Applied part | | Manufacturer |
| ⚠ | Caution, consult accompanying documents. | | |

SPECIFICATIONS

Type: Digital thermometer (Not predictive).

Range: 32,0°C - 42,9°C (90,0°F - 109,9°F) (°C / °F)

Accuracy: +/- 0,1°C, 35,5°C - 42,0°C (+/- 0,2 °F, 95,9 °F - 107,6°F).

Display: Liquid crystal display, 3 1/2 digits.

Memory: For storing the last measured value.

Battery: One 1.5 V DC, button battery (Size LR41 or SR41, UCC 392)

Battery life: Approx 200 hours

Dimensions: 12,4cm x 1,8cm x 0,9cm (L x W x H)

Weight: Approx 10 grams including battery

Using environment: Temperature: 5°C - 40°C (41°F - 104°F) Relative humidity: 15% - 95%RH

Storage and transportation condition: Temperature: -20°C - 55°C (-4°F-131°F)

Relative humidity: 15% - 95%RH

Ingress protection rating: IP 22

Classification: Type BF

DIRECTIONS

1. Press the On/Off button. A tone will sound as the screen will be lit, followed by the last recorded temperature. After showing the self-test temperature, the thermometer is now in the testing mode.

*Note: If the measured temperature is less than 32°C or 90°F, the LCD will display "Lo".
If the measured temperature is higher than 42,9°C or 109,9°F the LCD will display "Hi".

2. Position thermometer in desired location (mouth, rectum or armpit).
3. The degree sign flashes throughout the testing process. When flashing stops an alarm will beep for approximately 10 seconds. The measured reading will appear on the LCD simultaneously.
The minimum measurement time until the signaling tone (beep) must be maintained without exception. The measurements continues even after the buzzer notification. The displayed temperature will not change when thermometer is removed from its testing position.

*Note: Normally the buzzes are "Bi-Bi-Bi-Bi-"; Alarm beeps more rapidly when temperature reaches 37,8°C (100°F) or higher, and the buzzes are "Bi-Bi-Bi---Bi-Bi-Bi---Bi-Bi-Bi".

4. To prolong the battery life, press the On/Off button to turn off the device after testing complete.
If no action is taken, the unit will automatically shut off after around 10 minutes.
5. Store the thermometer in it's protective case.

BATTERY REPLACEMENT

1. Replace battery when "low battery" occurs in the lower corner of LCD-display.
2. Pull battery cover off as shown in Figure 1.
3. Gently pull out plastic circuit board with battery chamber approximately 1cm (See figure 2)
4. Use pointed object such as a pen to remove old battery. Discard battery lawfully.
Replace with new 1.5V DC button type LR41 or SR41, UCC392 or equivalent.
Be sure battery is installed with "+" polarity facing up. (See figure 3).
5. Slide battery chamber back into place and attach cover.

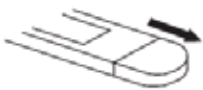
Figure 1

Figure 2

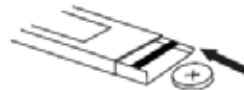
Figure 3

CLEANING AND DISINFECTION

Wipe the thermometer with a soft clean cloth. For stubborn stains, wipe the thermometer with a cloth that has been dampened with water or a neutral detergent solution and then wring thoroughly.
Finish by wiping with a soft dry cloth. Observe the following to prevent damage to the thermometer.

- △ Do not use benzene, thinner, gasoline or other strong solvents to clean the thermometer.
- △ Do not soak the sensing section in alcohol for long periods of time or attempt to sterilize it using hot water (Water at a temperature of 50°C (122°F) or higher).
- △ Do not use ultrasonic washing to clean the thermometer.

Disposal of this product and used batteries should be carried out in accordance with the national regulations for the disposal of electronic products.

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